

HPLC Application

ID No.: 20509

Metoprolol on Lux 5 μ m Amylose-2 in NP

Column: N/A

Dimensions: N/A

Order No: N/A

Elution Type: N/A

Eluent A: Hexane / EtOH (80:20) DEA (0.1%)

Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

Analyst Note: SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8470

- Holder Part No.: KJ0-4282

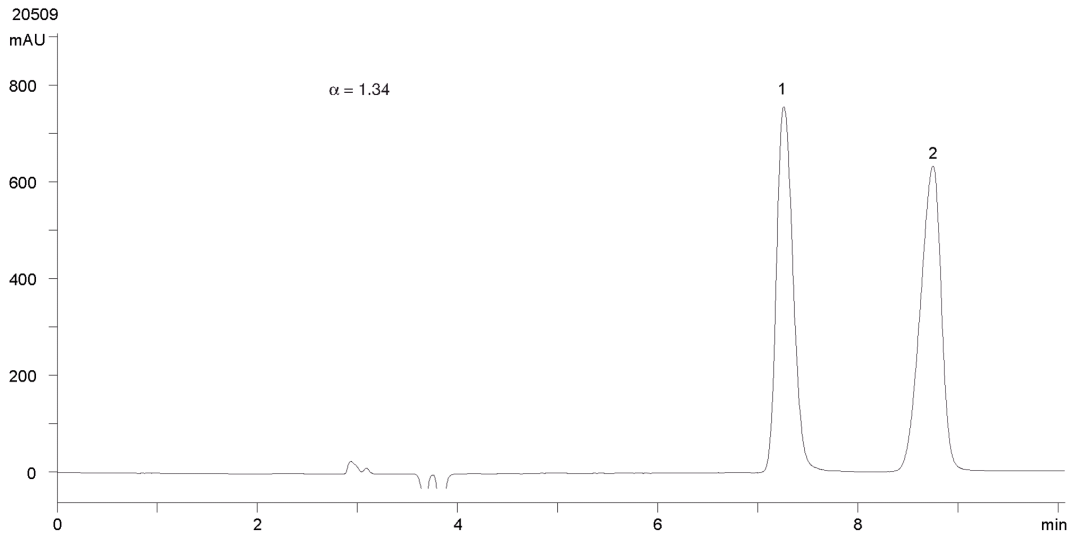

ANALYTES:

1 Metoprolol

Retention Time: 7.26 min

2 Metoprolol

Retention Time: 8.75 min

