

# HPLC Application

ID No.: 20512

## Tolperisone on Lux 5 $\mu$ m Amylose-2 in NP

**Column:** N/A

**Dimensions:** N/A

**Order No:** N/A

**Elution Type:** N/A

**Eluent A:** Hexane/ EtOH (98:2) DEA (0.1%)

**Gradient**

Step No.	Time (min)	Pct A
1	0	100

**Profile:**

**Flow Rate:** 1 mL/min


**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

**Analyst Note:** SecurityGuard<sup>™</sup> Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8470

- Holder Part No.: KJ0-4282

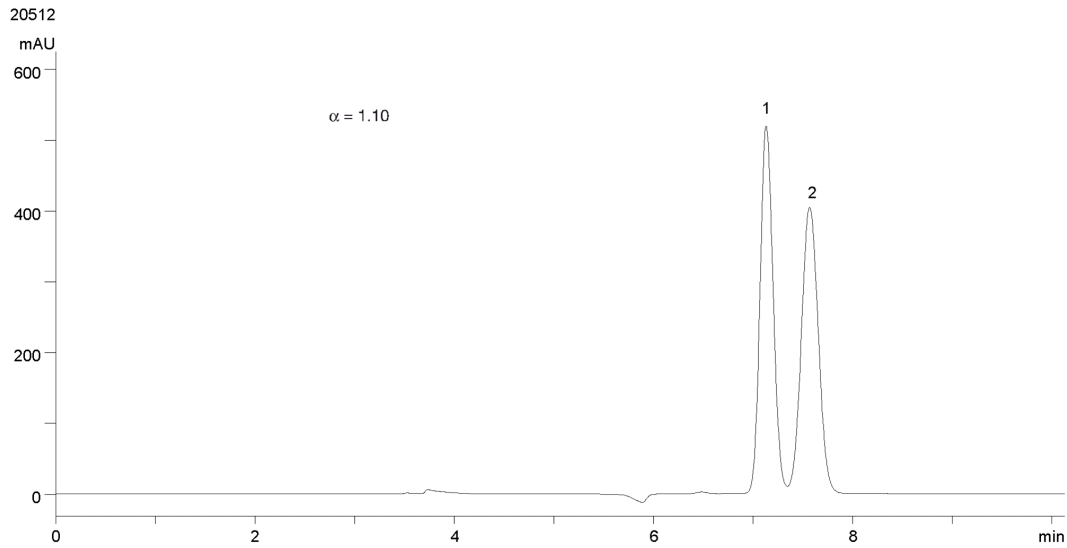

### ANALYTES:

**1** Tolperisone

Retention Time: 7.13 min

**2** Tolperisone

Retention Time: 7.57 min

