


Overlay: Yarra 3µm SEC-3000 versus TSKgel 5u G3000 SWxl high salt buffer MMQC Standard


Column: Yarra™ 3 µm SEC-3000, LC Column 300 x 7.8 mm, Ea
Dimensions: 300 x 7.8 mm ID
Order No: 00H-4513-K0
Elution Type: Isocratic
Eluent A: 50 mM Sodium Phosphate pH 6.8 / 0.3 M Sodium Chloride
Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

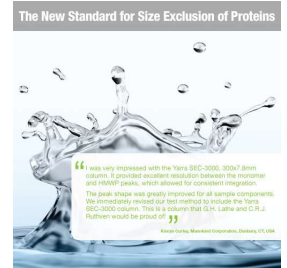

Products used in this application:

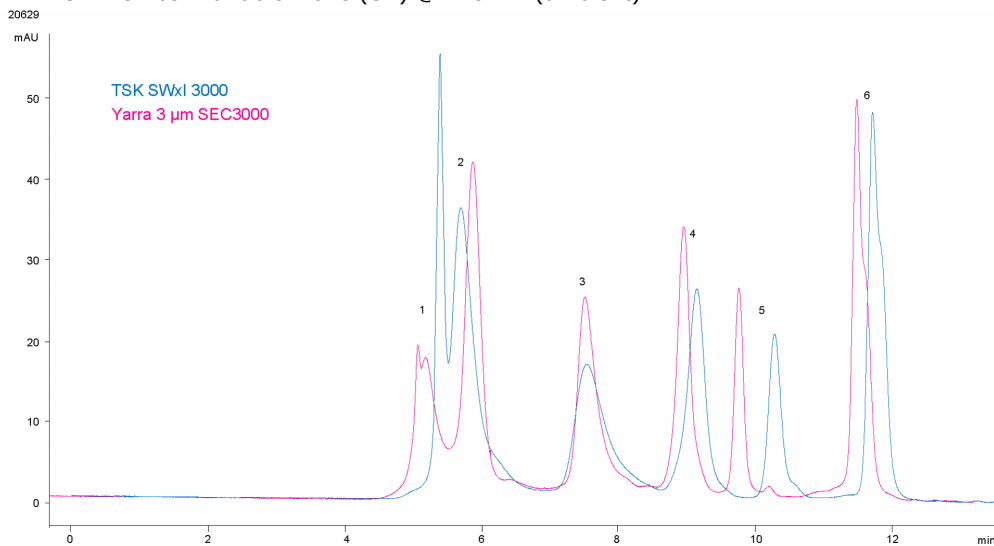

ANALYTES:

- 1** IgM (900kDa)
Retention Time: 5.5 min
- 2** Thyroglobulin (669kDa)
Retention Time: 6.3 min
- 3** Beta Amylase
Retention Time: 7.9 min
- 4** Ovalbumin (44kDa)
Retention Time: 9.3 min
- 5** Myoglobin (17kDa)
Retention Time: 10.1 min
- 6** Uridine
Retention Time: 11.7 min

