

HPLC Application

ID No.: 20702

Lansoprazole on Lux 5 μ Cellulose-1 in NP Elution Mode (0.1% DEA)

Column: Lux[®] 5 μ m Cellulose-1, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4459-E0

Elution Type: Isocratic

Eluent A: 0.1% DEA IN HEXANE

Eluent B: 0.1% DEA IN IPA

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	20	60	40

Flow Rate: 1


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

20702



Products used in this application:

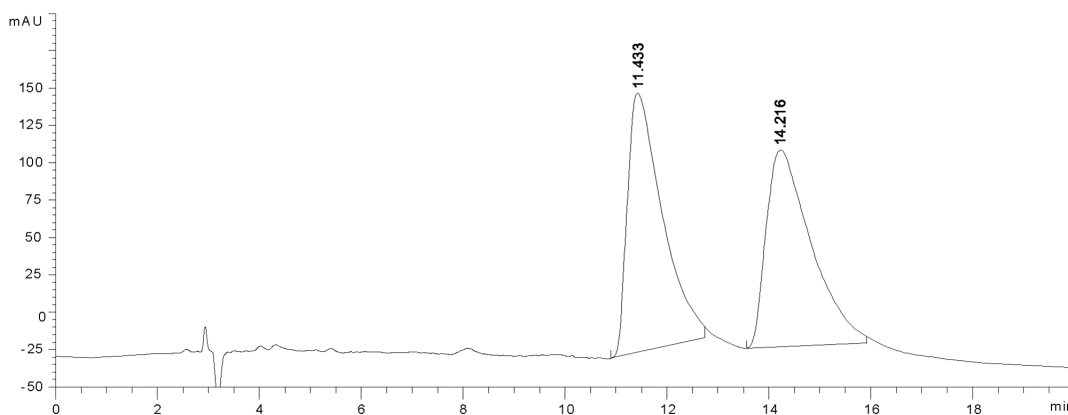

ANALYTES:

- 1 Lansoprazole
- 2 Lansoprazole

