

HPLC Application

ID No.: 20711


Rabeprazole on Lux 5 μ m Amylose-2 in PO Elution Mode (0.1% FA)

Column: N/A
Dimensions: N/A
Order No: N/A
Elution Type: N/A
Eluent A: 0.1% FA in CH₃CN
Eluent B: 0.1% FA in IPA
Gradient

Step No.	Time (min)	Pct A	Pct B
Profile: 1	20	95	5

Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

19711

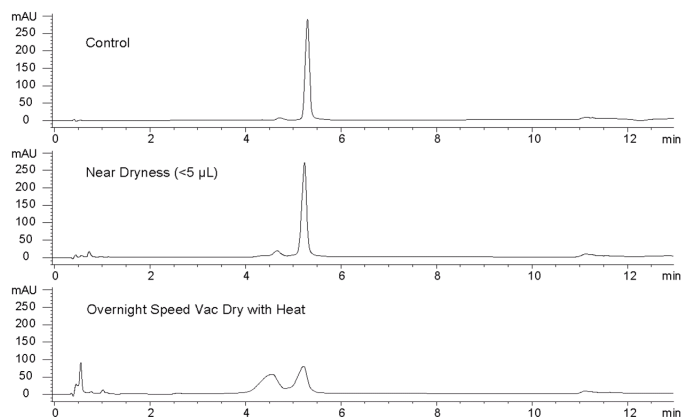

ANALYTES:

- 1 Rabeprazole
- 2 Rabeprazole

