

HPLC Application

ID No.: 21806

Loctofen on Lux 5 μ m Cellulose-2 in RP


Column: Lux[®] 5 μ m Cellulose-2, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4457-E0
Elution Type: Isocratic
Eluent A: MeOH / 20mM NH₄HCO₃+0.1%DEA (90:10) with DEA (0.1%) in buffer
Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)



Products used in this application:

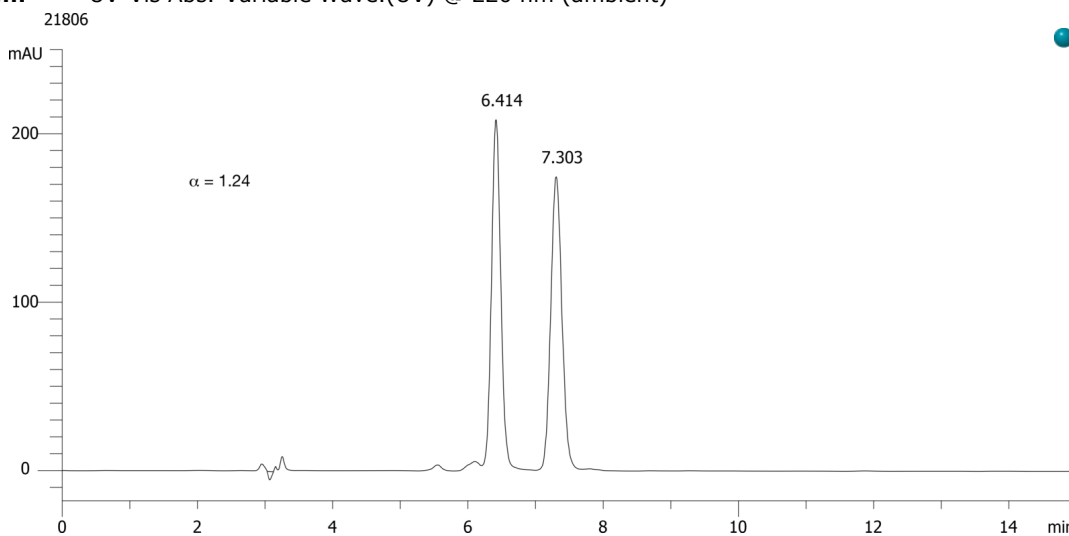

ANALYTES:

- 1 Loctofen
Retention Time: 6.414 min
- 2 Loctofen
Retention Time: 7.303 min