

HPLC Application

ID No.: 22724

Sufentanil at high pH using Kinetex 5u EVO C18 50 x 2.1 mm

Column: Kinetex® 5µm EVO C18 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4633-AN

Elution Type: Gradient

Eluent A: 10mM Ammonium bicarbonate

Eluent B: Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	0.5	90	10
	3	2	75	25
	4	4.5	20	80
	5	4.51	15	85
	6	5.5	15	85
	7	5.51	90	10
	8	7	90	10

Flow Rate: 0.5 mL/min

Col. Temp.: 40 °C

Detection: Tandem Mass Spec (MS-MS) @ (ambient)


Detector Info:

22724

SCIE<



Products used in this application:

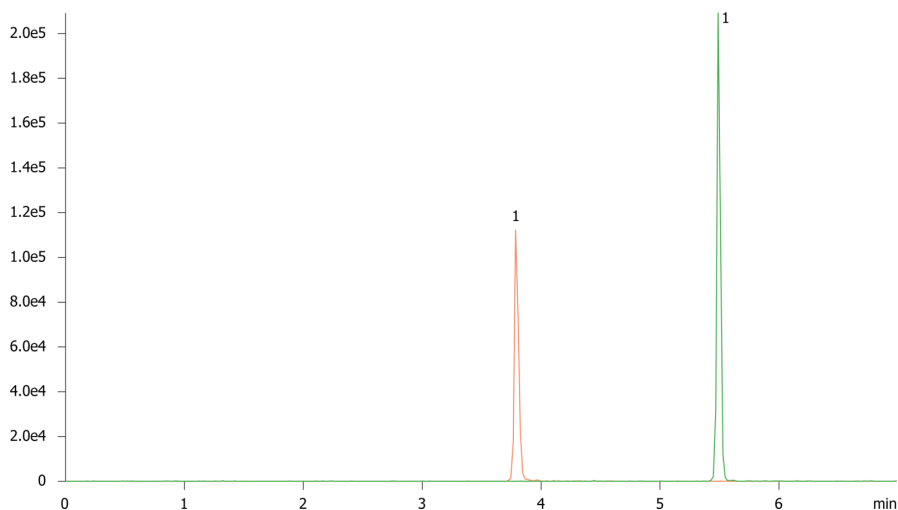

ANALYTES:

1 Sufentanil

Retention Time: 3.79 min

