

HPLC Application

ID No.: 22818

Trace API Impurities - Kinetex 1.7u F5 50x2.1mm

Column: Kinetex[®] F5 1.7 μ m 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4722-AN

Elution Type: Gradient

Eluent A: 20mM Potassium phosphate pH 2.3

Eluent B: Methanol


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	10	5	95
	3	10.1	95	5

Flow Rate: 0.3 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

22818

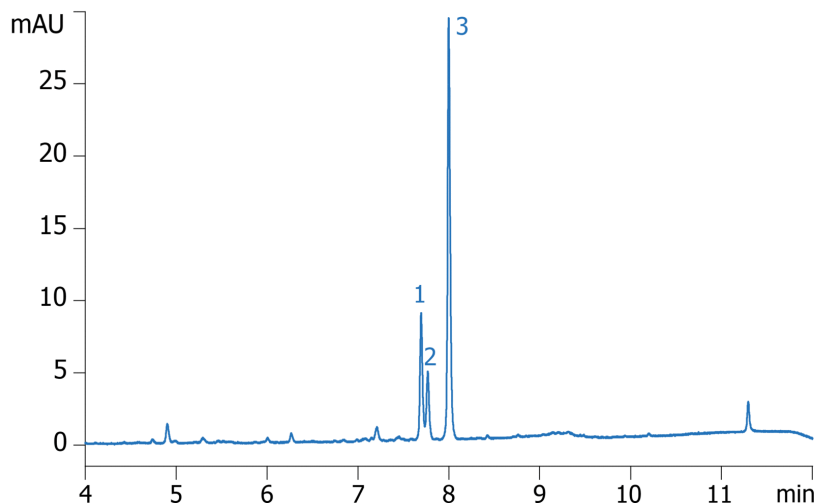

ANALYTES:

- 1 Impurity 1
- 2 Impurity 2
- 3 API



Products used in this application:

