

HPLC Application

ID No.: 23753

Zopiclone on Lux 5µm i-Cellulose-5 250x4.6mm (UV)

Column: Lux[®] 5 µm i-Cellulose-5, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4756-E0

Elution Type: Isocratic

Eluent A: 5mM Ammonium Bicarbonate

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	20	80


Flow Rate: 1 mL/min

Col. Temp.: 20 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 300 nm (20 °C)
23753



Products used in this application:

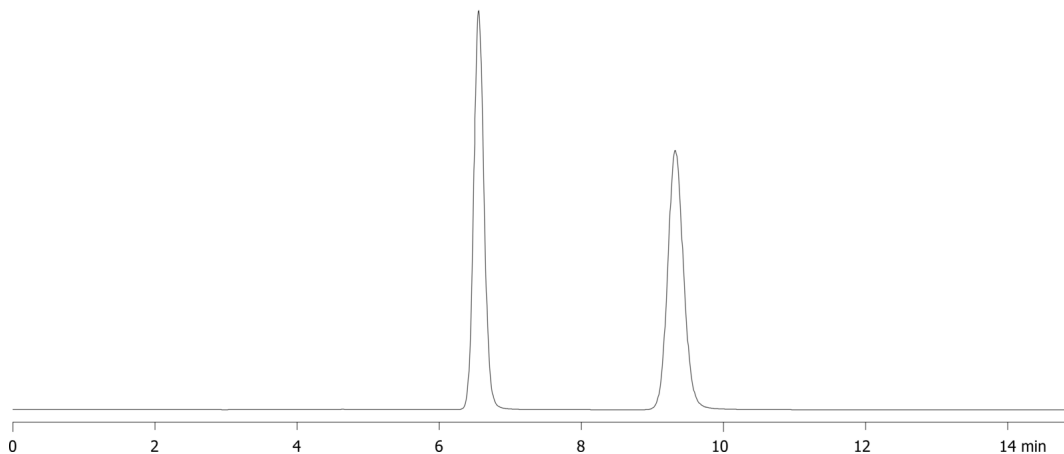

ANALYTES:

1 Zopiclone

Retention Time: 6.552 min

2 Zopiclone

Retention Time: 9.319 min

