


## Rituximab on bioZen 1.8 µm SEC-3

Column: bioZen<sup>™</sup> 1.8 µm SEC-3 , LC Column 300 x 4.6 mm, Ea  
Dimensions: 300 x 4.6 mm ID  
Order No: 00H-4772-E0  
Elution Type: Isocratic  
Eluent A: 50 mM K<sub>2</sub>PO<sub>4</sub> 100 mM Na<sub>2</sub>SO<sub>4</sub>, pH 6.8  
Gradient  
Profile: 

Step No.	Time (min)	Pct A
1	30	100

  
Flow Rate: 0.2 mL/min  
Col. Temp.: 25 °C  
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (25 °C)

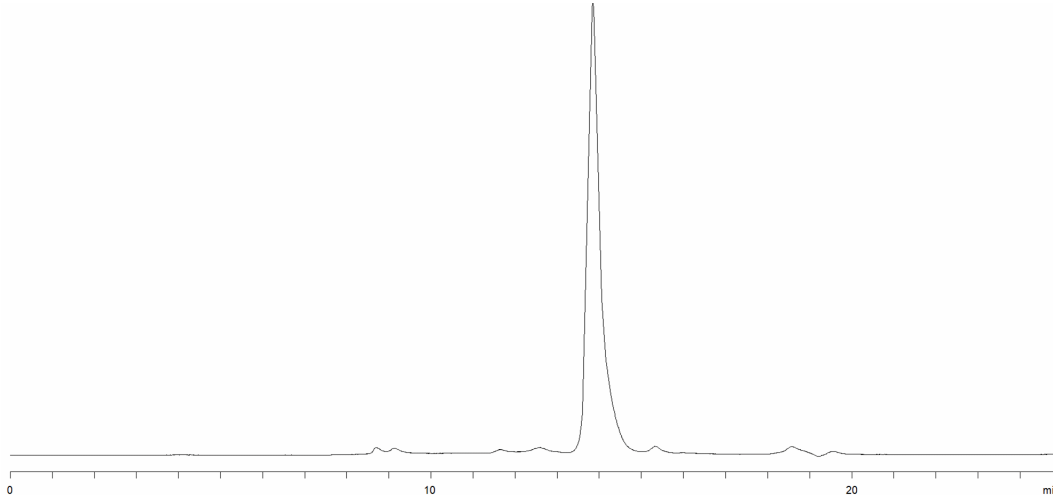

### ANALYTES:

1 Rituximab

