

HPLC Application

ID No.: 25042

Trastuzumab on bioZen 1.8 μ m SEC-3

Column: bioZen[™] 1.8 μ m SEC-3 , LC Column 300 x 4.6 mm, Ea

Dimensions: 300 x 4.6 mm ID

Order No: 00H-4772-E0

Elution Type: Isocratic


Eluent A: 50 mM KH₂PO₄ 250 mM KCl, pH 6.8

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 0.35 mL/min

Col. Temp.: 30 °C

Detection: LC/UV (DAD, PDA) @ 280.000000000 nm (nanometers) (ambient)

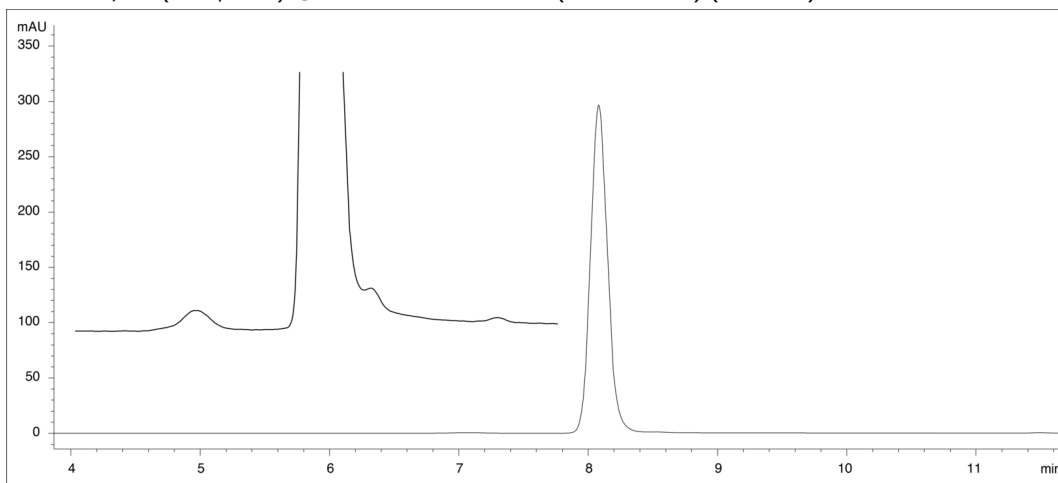

ANALYTES:

- 1 Trastuzumab
- 2

