

HPLC Application

ID No.: 25060

Charge Variants of Cetuximab Biosimilar using pH Gradient on bioZen 6 µm WCX

Column: bioZen[™] 6µm WCX, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4777-E0

Elution Type: Gradient

Eluent A: CX-1, pH 5.6


Eluent B: CX-1, pH 10.2

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	100	0
	2	1	100	0
	3	21	0	100
	4	23	0	100
	7	24	100	0

Flow Rate: 1 mL/min

Col. Temp.: 30 °C

Detection: LC/UV (DAD, PDA) @ 280.000000000 nm (nanometers) (ambient)

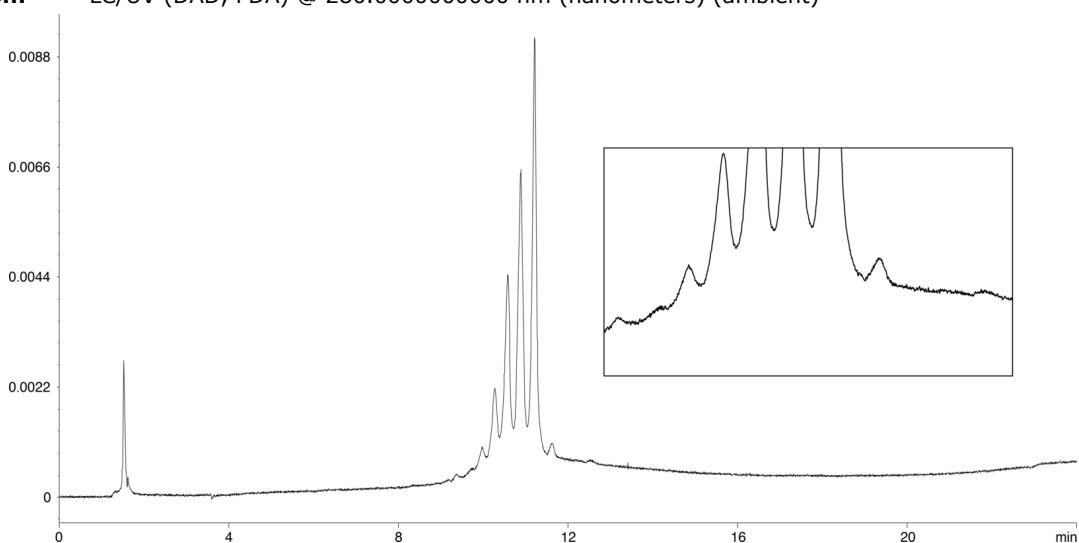

ANALYTES:

1 Cetuximab

